

PRELIMINARY:



EXPERIMENT	OBSERVATION	INFERENCE
Colour of the salt is noted	white/colorless	Absence of $\text{Cu}^{2+}/\text{Ni}^{2+}/\text{Mn}^{2+}$
Odour is noted	a) vinegar like smell	Presence of acetate (CH_3COO^-)
	b) smell of ammonia	Presence of NH_4^+ ion
	c) No characteristic odour	Absence of NH_4^+ or CH_3COO^-
Solubility is noted	1) Soluble in water	Presence of NH_4^+ and or Cl^-
	2) Insoluble in water	Absence of NH_4^+ or Cl^-

TEST FOR Cl^- , Br^- , I^-

Take a dry test tube, add a pinch of salt into it. Then add 1-2 ml of conc. H_2SO_4 to it and heat.

a) Colourless gas which gives a dense white fumes when a glass rod dipped in NH_4OH is shown near the mouth of test tube

Presence of Cl^-

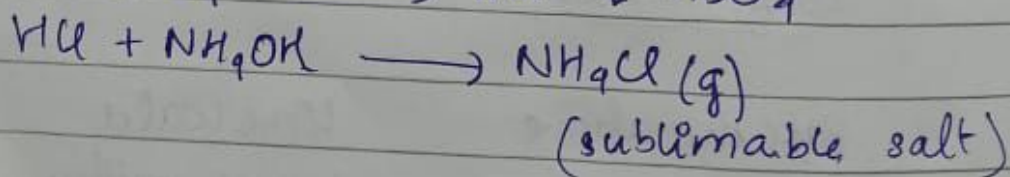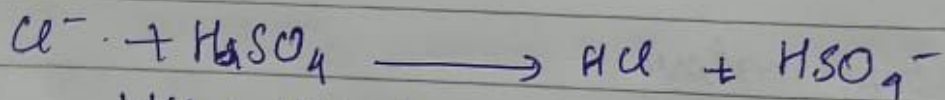

Confirmatory Test

To the salt solution add ^{1-2 ml} dil. HNO_3 . Then heat it and then cool it. Now add AgNO_3 to it

white curdy precipitate which dissolves in NH_4OH soln.

Presence of Cl^- is confirmed


TEST FOR ACID RADICAL (Anion)

Test for CO_3^{2-} (Carbonate ion)

Take aml of
dil. HCl in test
tube then add
a pinch of
salt into it.

Brisk
effervescence
with the
liberation
of colourless,
odourless
gas



Presence of
 CO_3^{2-}

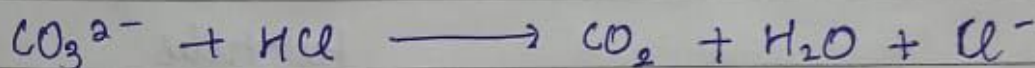

Confirmatory:


Pass the above
gas through
lime water

lime water
turns milky

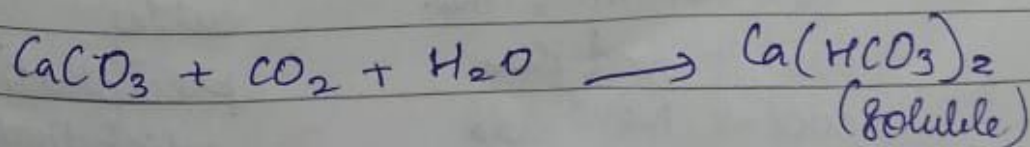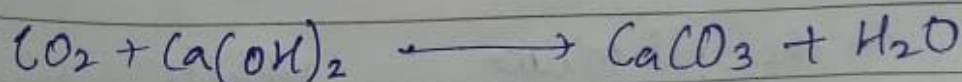


Cl⁻ - greenish
I⁻ - violet
P⁻ - pale yellow

Page No.:
Date:

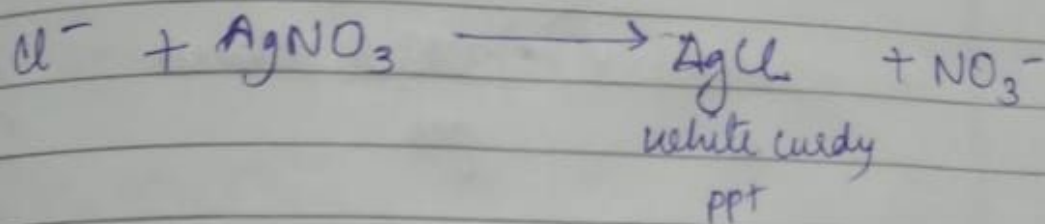

b) Reddish brown gas Presence of Br⁻

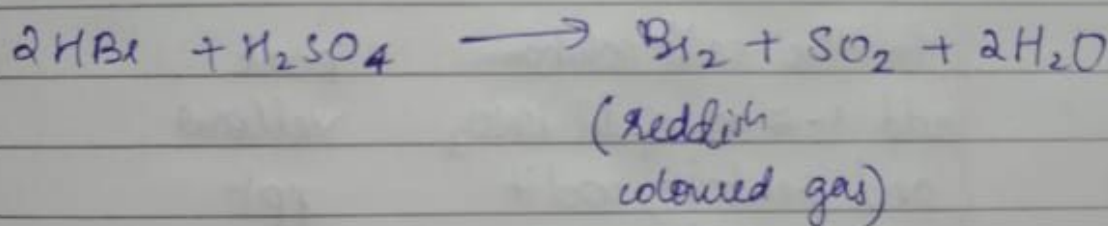


Confirmatory Test:

To the salt soln
add ~~the~~ dil HNO₃,
then heat and
cool it, Then add
AgNO₃ to it

Pale yellow ppt



Presence of Br⁻ is confirmed

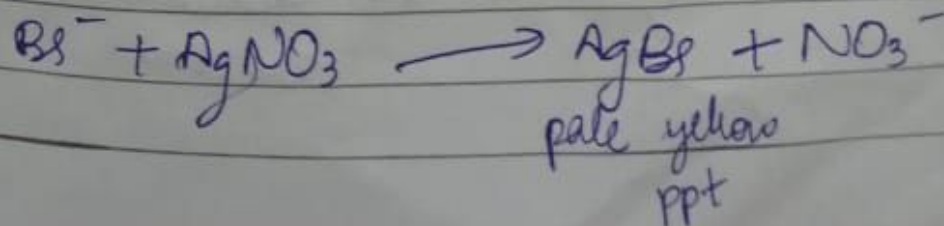

X only for viva

c) Violet coloured gas

Presence of I^- .

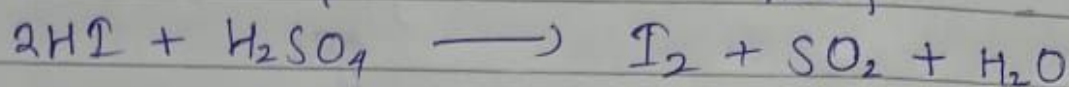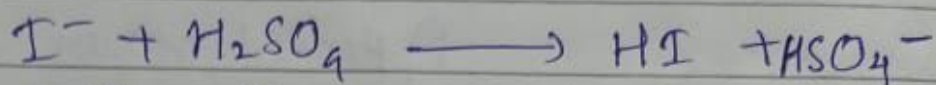


Confirmatory:

To the salt solution
add 1-2 ml of dil. HNO_3
and heat it, cool it
Then add $AgNO_3$.

dark
yellow
ppt

presence of
 I^- is
confirmed

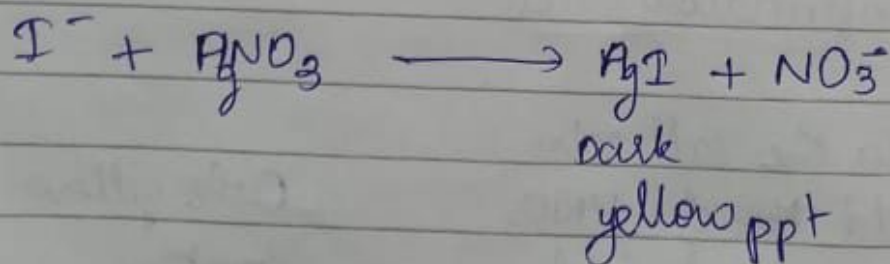

or

No characteristic
reaction



Absence
of Cl^- ,
 Br^- , I^-

TEST FOR NITRATE (NO_3^-)

To the above mixture add paper balls and heat strongly

Plenty of pale brown gas is produced

Presence of NO_3^-

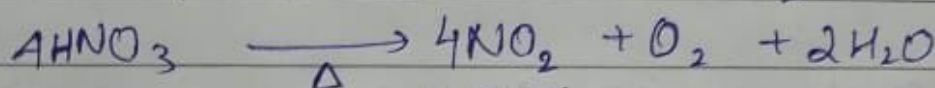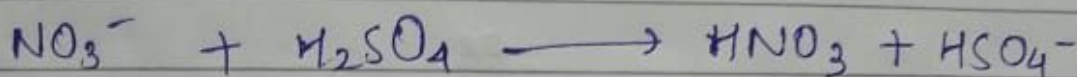

(pale brown coloured gas)

Confirmatory Test (Brown ring test)

Take salt soln in the test tube (saturated)
Add equal volume of freshly prepared FeSO_4 soln.
mix it well, then add conc. H_2SO_4 along the sides of the test tube without shaking


Brown ring is formed

Presence of NO_3^- is confirmed.

M

Page No.:

Date:

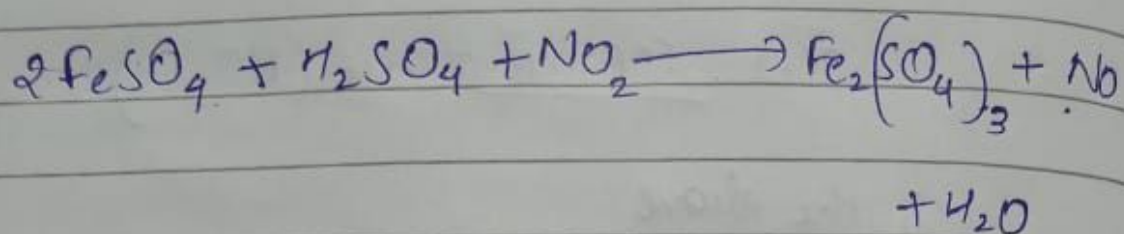

nitroso ferrous sulphate
(Brown ring)

or

Dark

No characteristic
reaction



Absence of
 NO_3^-

TEST FOR SULPHATE ACETATE

Take salt
in a watch
glass add
equal quantity
of oxalic acid,
then add 1-2
drop of water
Rub it using
your finger
and smell it

Vinegar like
smell

Presence
of acetate
 CH_3COO^-

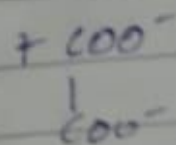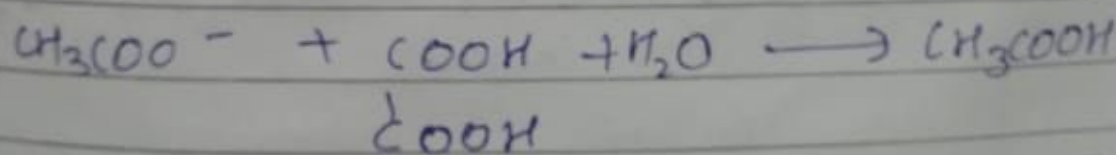

Confirmatory

To the salt
soln add
neutral FeCl_3 soln

Red blood
colouration

Add NH_4OH drop by drop to $\text{FeCl}_3 \Rightarrow$ neutral

Page No: FeCl_3 .

Date:

The red blood coloured soln is divided into two parts.

To one part add dil HCl .



Red colour disappears.

Presence of CH_3COO^- confirmed

To the second part add water and boil it

Reddish brown ppt

Presence of CH_3COO^- confirmed.

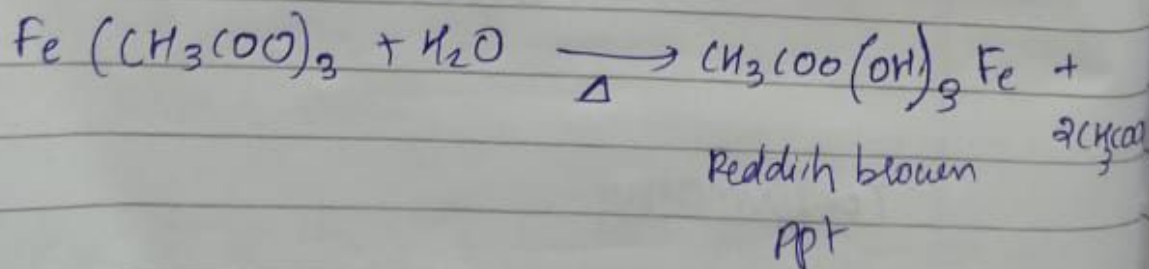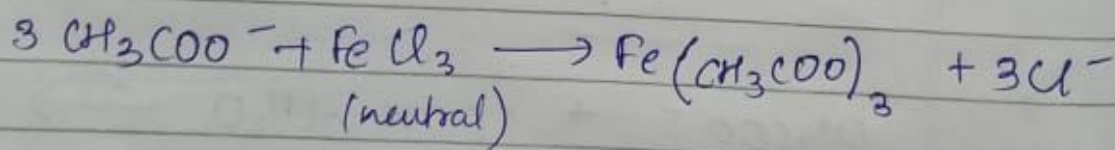


Test for SO₄ SULPHATE (SO_4^{2-})

To the salt
soln add
 BaCl_2 soln.

white ppt
insoluble in
conc. HCl .

Presence of
 SO_4^{2-}

Presence
of CH_3COO^-
confirmed

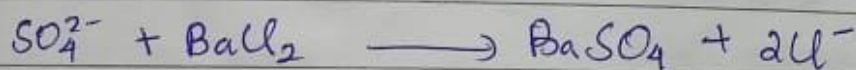

Presence
of CH_3COO^-
confirmed-

Confirmatory


To the salt
soln add lead
acetate soln

white ppt
soluble in
excess of
ammonium
acetate soln

presence
of SO_4^{2-}
confirmed

+ 3Cl^-

$\text{Pb} + 2\text{CH}_3\text{COO}^-$
when

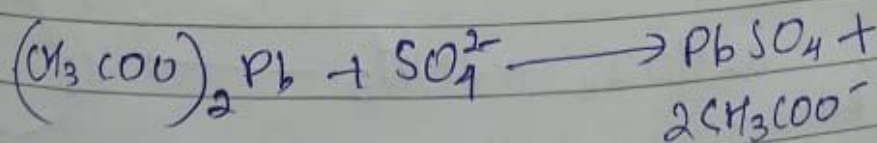

TEST FOR BASIC RADICAL


Test for O group:

To the dry salt
add conc. NaOH
and heat

Smell of
ammonia

Presence of
 O group
(NH_4^+)

confirmatory:

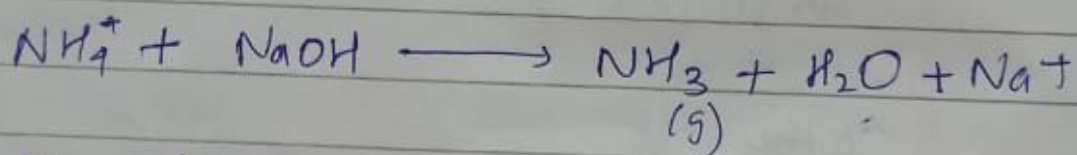


confirmatory:




Pass the above
gas through
Nessler's reagent

Brown ppt

Presence of
 NH_4^+ Confir.

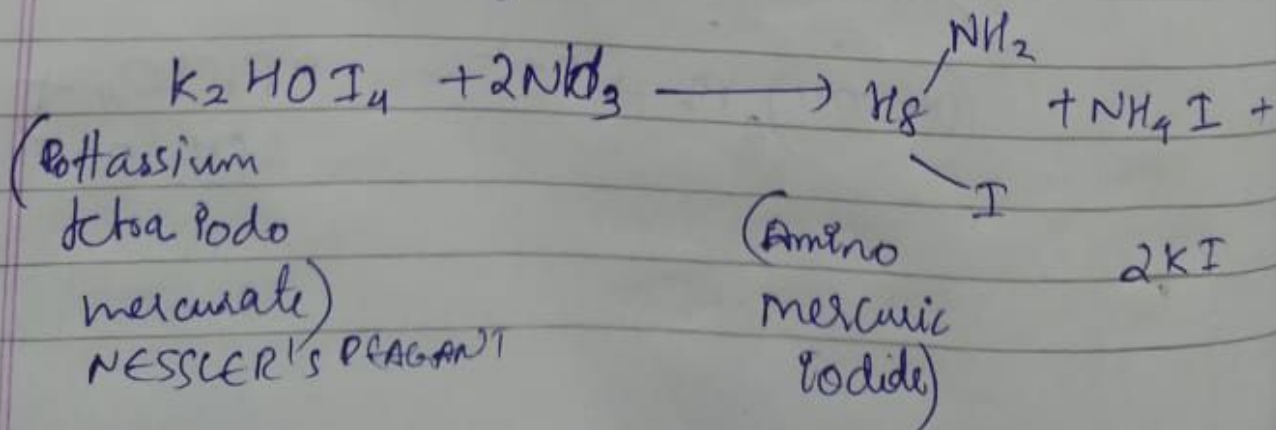

Group reagent

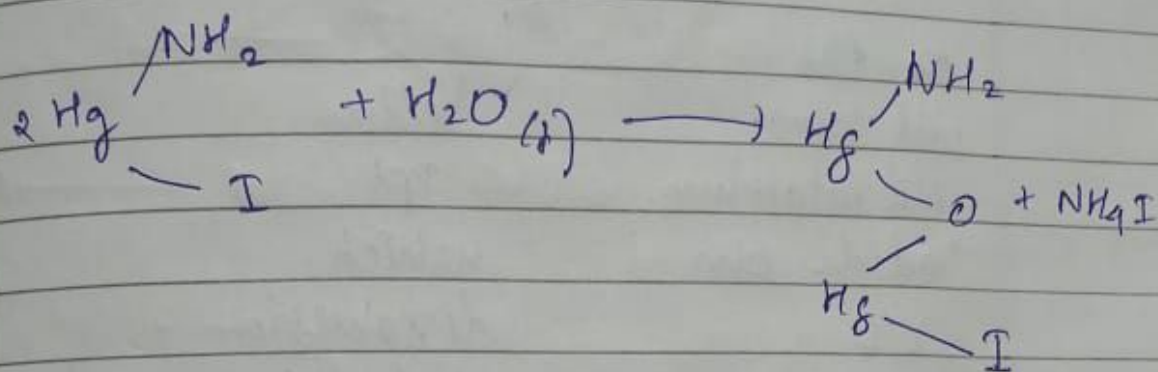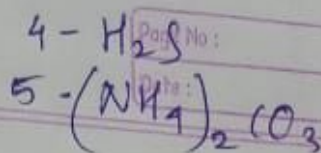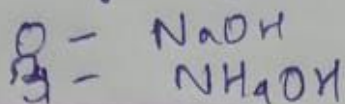


Test for 1 group:

To the salt soln add dil HCl

white ppt

presence of group 1.

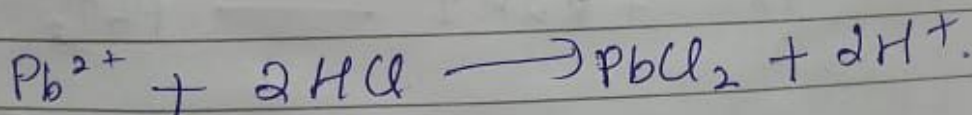


Confirmatory:

To the salt soln add Potassium chromate K_2CrO_4

beige yellow ppt


Presence of Pb^{2+} confirmed

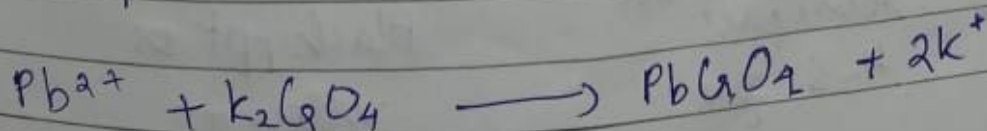

To the
salt soln.
add potassium
iodide soln
KI


Yellow
ppt
which
dissolves in
hot water
and reappears
as golden
sparkles
on cooling

Presence
of Pb^{2+}
confirmed

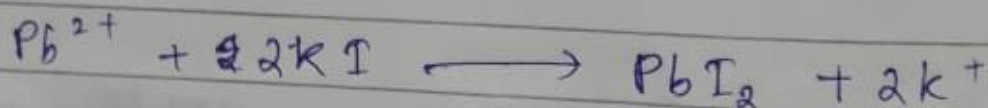

X only for viva

Test for group 2

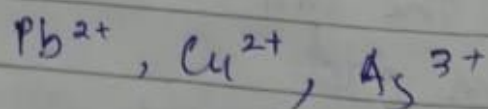

To the salt
soln, pass
dil HCl
 H_2S gas

Black ppt or
Yellow ppt

Presence of
group 2.


Test for group 3 Al^{3+} , Fe^{3+} X viva

To the salt

soln add excess
of solid NH_4Cl
then add
excess of NH_4OH

gelatinous
white ppt

Presence of
group 3
(Al^{3+})



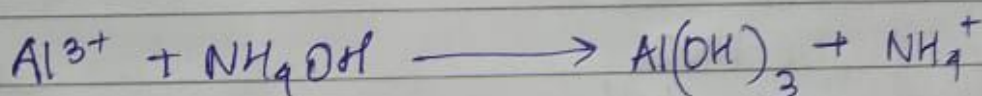
Confirmatory:

(Blue lake test)

Dissolve the
above gelatinous
ppt in dil. HCl .
Add blue litmus
soln. Then add
 NH_4OH soln
in excess
drop by drop

Blue lake
formation

Presence of
 Al^{3+}
confirmed

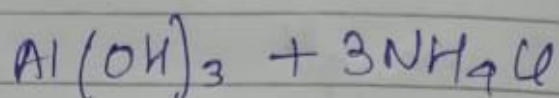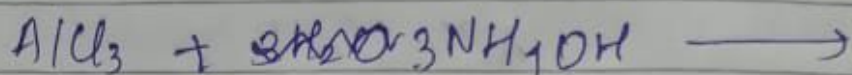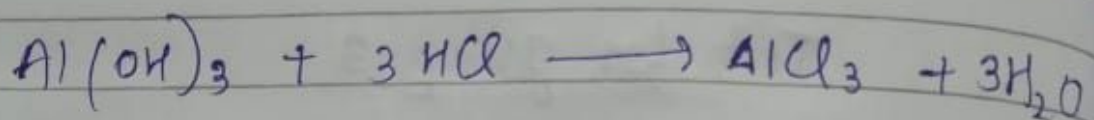


Blue ash test

to the salt
soln add
dil HNO_3
and cobalt
nitrate soln
then heat
it, then dip a
filter paper into
the soln, burn
the filter paper

Blue
ash is
formed

Presence of
 Al^{3+}
confirmed

Test for 4 group

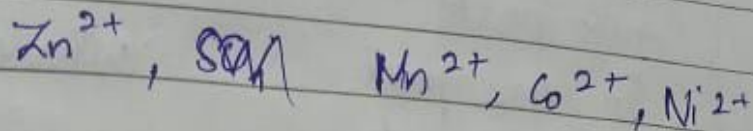


To the salt
soln add solid
 NH_4Cl and
add excess
of NH_4OH
then add
 H_2S

dirty
white ppt

Presence
of group 4
(Zn^{2+})

confirmatory

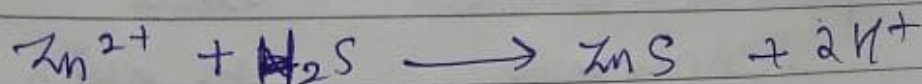

Confirmatory :



To the
salt so
dissolve the
white ppt
in dil HCl
boil off H_2S
then chloride

the soln into
two parts.

To one part
add NaOH soln
dropwise

white ppt


Presence
of Zn^{2+}
confirmed

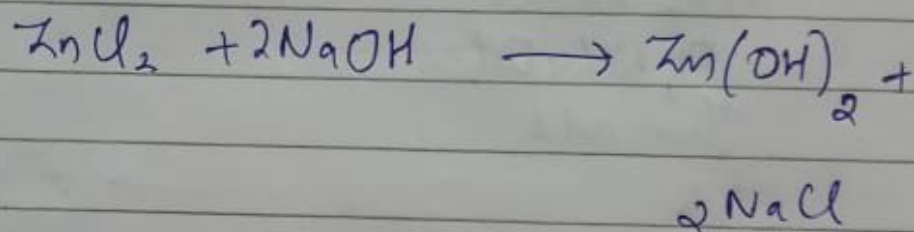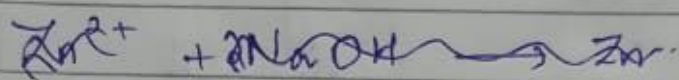

To the second
part add Potassium
ferrocyanide


white
or bluish
white ppt

Presence
of Zn^{2+}
Confirmed

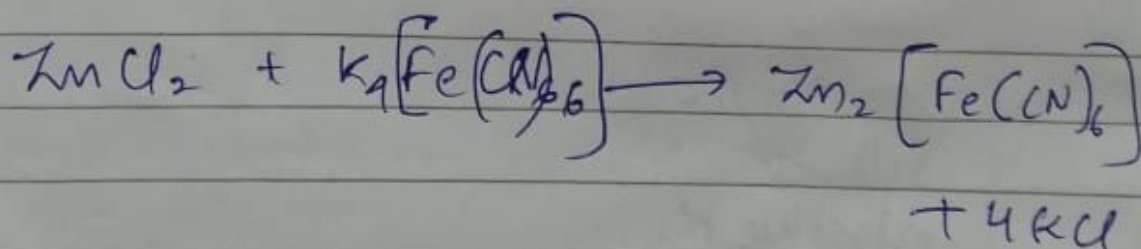

Page No :
Date :

Test for group 5:



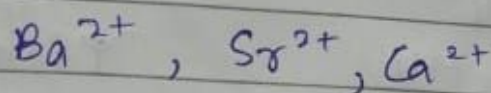


To the salt
soln add
solid NH_4Cl
add excess

white ppt

presence
of group 5

of NH_4OH
then add
ammonium
carbonate soln

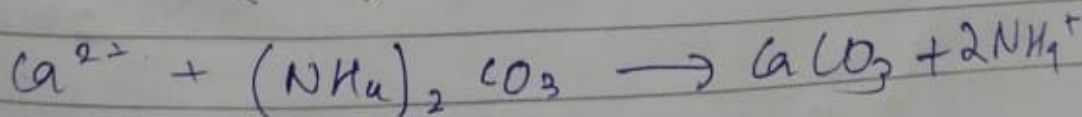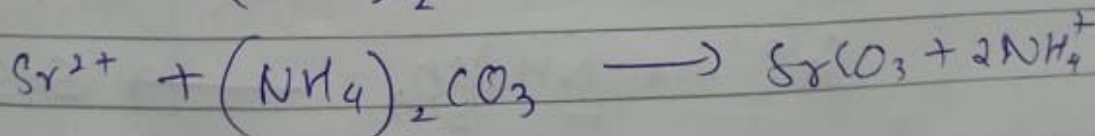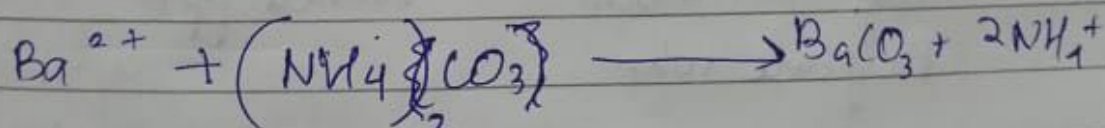


confirmatory

redissolve the
above ppt in
acetic acid
then boil off CO_2

Then divide into
3 parts
To the
first part add
Potassium Chromate
(K_2CrO_4)

Yellow
ppt


Presence of
 Ba^{2+}
confirm

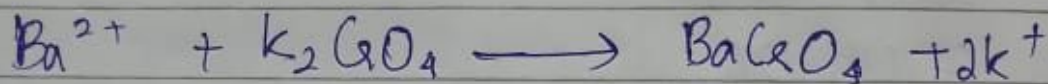

To the
second part
add ammonium
sulphate soln

White ppt

Presence of
 Sr^{2+}
Confirmed

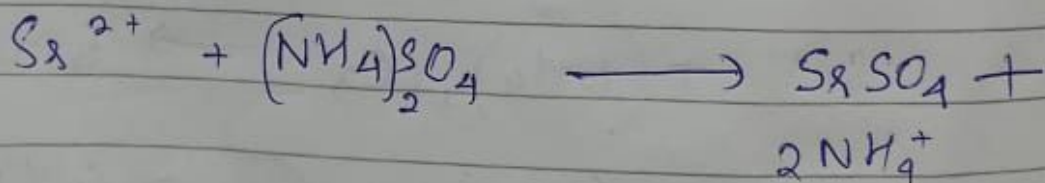


To the third
part add ammonium
oxalate soln.

White ppt
while scratching
the side
of the test
tube

Presence
of Ca^{2+}
confirmed

conc HCl is used to convert it to chloride bcz it is more volatile.

Date:

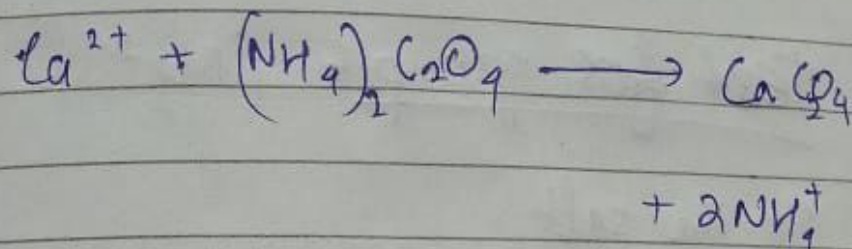

Flame test:

~~et~~ e^- transition is ^{from lower level to higher level and vice versa} the principle

Take a dry watch glass add pinch of salt into it. Then add

2-3 drops of ~~conc. HCl~~ ^{conc. HCl} and make

paste of salt. Take the salt in the glass rod. Show it to the flame.

Note the colour of flame.

1) apple green

Ba^{2+}

confirmed

2) crimson red

Sr^{2+}

confirmed

3) Brick red

Ca^{2+}


confirmed

Test for group 6 (Mg^{2+})

To the salt
soln add solid
 NH_4Cl and
then excess of
 NH_4OH . Then
add disodium
hydrogen phosphate
 Na_2HPO_4

white ppt

Presence
of group
6.

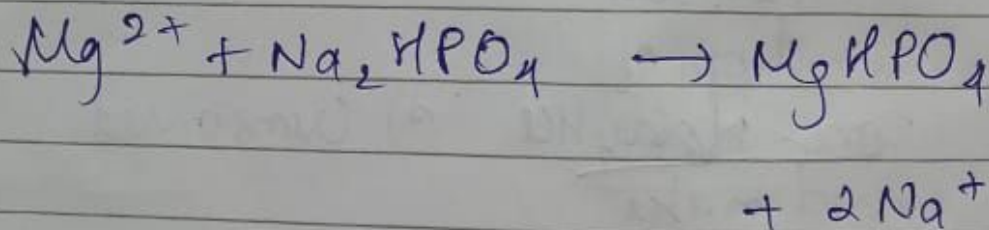

Confirmatory

To the salt
soln add magneson
reagent

Blueish
white
ppt

Presence
of Mg^{2+}
confirmed

Organic compounds

papergrid

Date: / /

1. Preliminary Test
Colour noted Colourless absence of phenol, cresol etc
2. Odour
Carbolic smell Presence of phenolic grp.
Pleasant fruity smell Presence of Ester

No characteristic odour Absence of Phenol/Ester

3. Take the given organic compound & burn it Sooty flame aromatic
non-sooty flame aliphatic

4. Test for unsaturation

Baeyer's test

To the org. compound soln. Violet colour ^{disappears} unsaturated
add alk. $KMnO_4$ Violet colour persists saturated


Test for Alcohol


Take the soln. of org. compound & dip a blue litmus paper. turn red presence of alc. group

blue litmus paper.

5. Reaction with Na

Add a piece of Na metal Evolution of gas Presence of alcoholic group

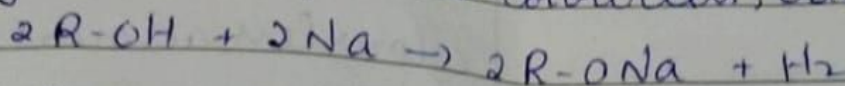

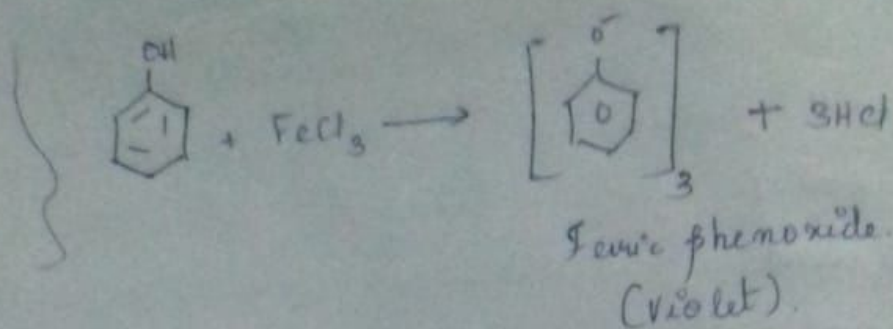


3) Test for carboxylic acid. $>\text{C}=\text{O}$
(Litmus test)

→ NaHCO_3 test

i) To the organic compound
solid sodium bicarbonate.

Brisk effe-
-escence with
the liberation of
colorless, odourless
gas

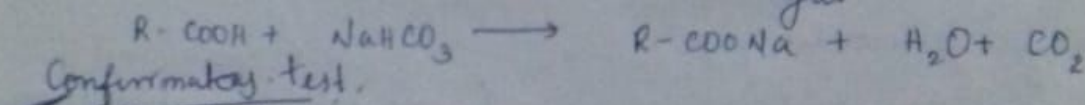

Confirmatory test.

ii) Ester Test

i) To the organic compound.
with ethanol. $\text{C}_2\text{H}_5\text{OH}$

4) Test for carbonyl grp.
2,4 DNP test

→ To the given
organic compound.
add 1 ml alcohol

Bayer's test - test for saturation/unsaturation.

• Take KMnO_4 sol.
make it alkaline

→ To the soln of
organic compound.
add alkaline KMnO_4 .

- Violet color disappears

unsat
- unsat.

- violet color persists

saturated

1) Test for alcohol. (Aliphatic saturated)


Take the soln of
organic compound,
dip in blue litmus
paper.

litmus turns
red.

Presence of
OH grp

→ Add a piece of Na metal
to organic compound.

Presence of
OH grp

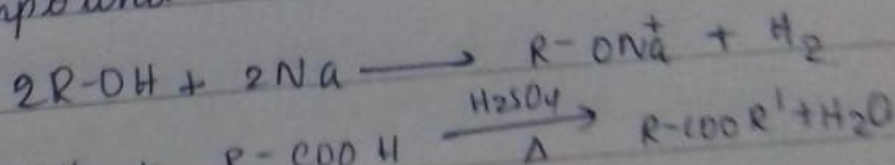

Confirmatory grp


→ Take organic compound,
then add 2ml of ethanol

acid into it
→ Add 1 or 2 ml of
conc H_2SO_4 and
heat strongly.

→ Transfer the contents
into 15-20 ml of 5%
aq. Na_2CO_3 soln. taken
in a beaker. Smell the
content.

pleasant fruity
smell.

presence
of OH
- grp.
confirmed

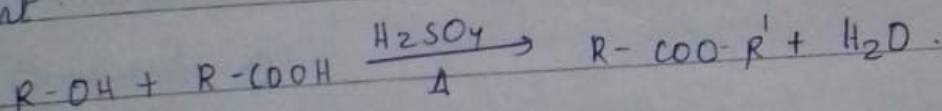

2) Test for phenolic grp.

→ Blue litmus test

turns red.

presence
of phenolic
grp Ar-OH.

1. To the soln. of organic
compound add blue litmus


To the aq. soln of organic
compound add neutral
ferric chloride soln.

violet coloration.

Presence
of phenol
grp.

then add 4 ml of
alcohol

2, 4 DNP reagent a yellow, orange
and shake the or red ppt
mixture. Will separate out

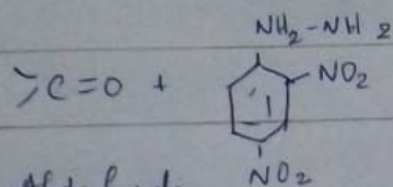

Test for Aldehyde

To the organic compound
add Schiff's base.

pink color

presence of
aldehyde

Confirmatory

- i) Tollen's reagent is
added to organic compound.
Heat strongly in a water bath.

appearance of
~~the~~ silver mirror
or grey ppt

Presence
of CHO
grp.

Take equal volume of
Fehling's reagent no. 1 & 2.
Then add organic compound
and heat strongly in
water bath.

Reddish brown
ppt.

Prese
nce of
-CHO.

Test for ketone

Barros Sodium Nitroprussiate test

Take 1 ml of Sodium nitroprusside soln. and make it alkaline by adding a few drops of NaOH soln. \rightarrow To it add a pinch of Organic compound and shake it

A red blood coloration

Presence of ketone

Meta dinitro benzene test

To the given organic compound add 1 ml of ethanol and then add 1 ml of ethanolic soln of Metadinitro benzene. and a few drops of NaOH soln. and then shake.

Violet coloration

presence of t-c-i